


Micro-Touch[®] Affinity[®]

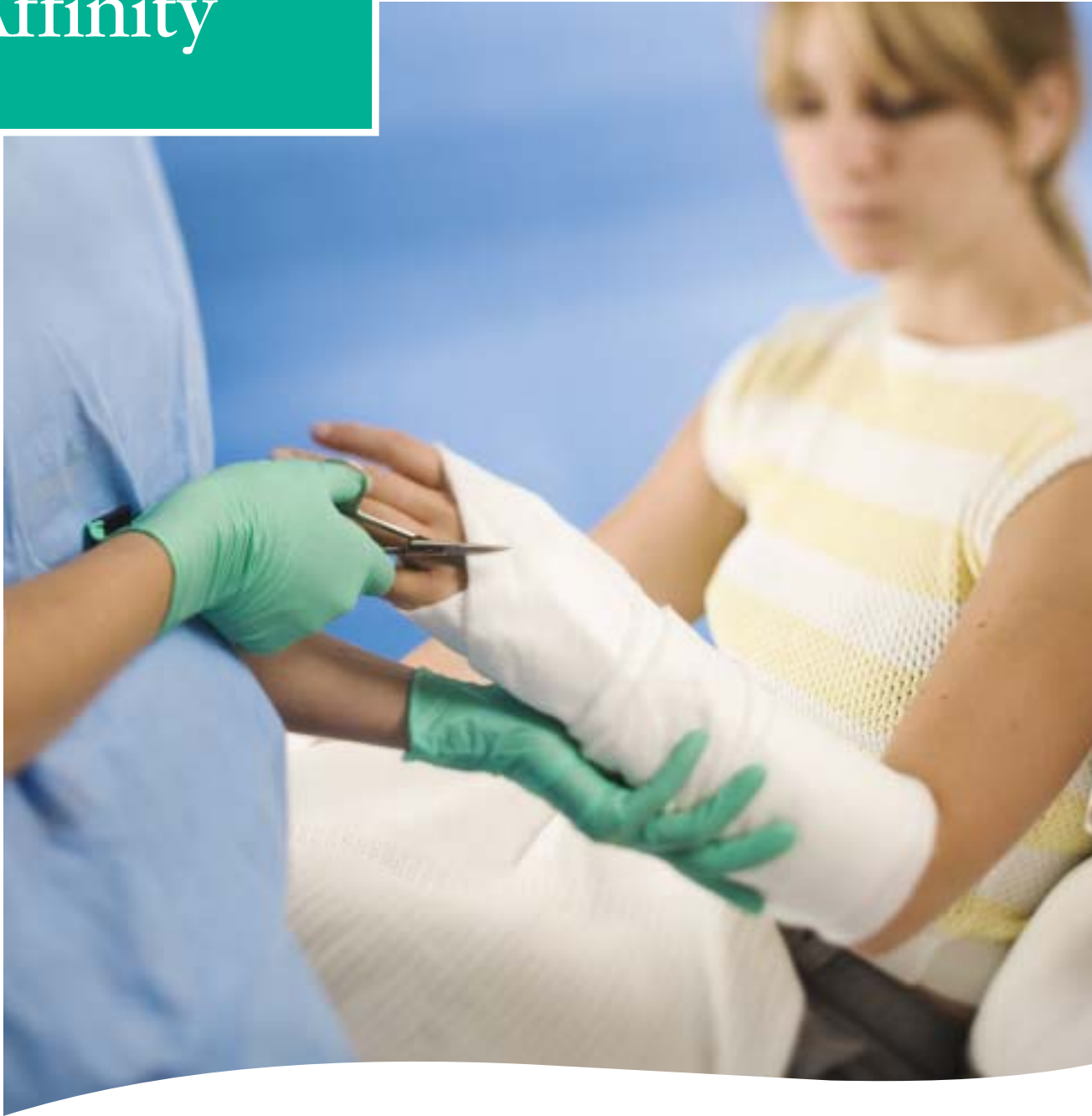

Neoprene Medical Examination Glove
True comfort and grip

Ansell

Micro-Touch® Affinity®

Latex Free and Powder Free Medical Examination Glove

At Ansell, we are always looking to take medical glove technology to the next level, and that is precisely what we have done with our Micro-Touch Affinity examination gloves. Created from neoprene, Micro-Touch Affinity offers the fit, comfort, protection, and stretchability you are looking for.

How does neoprene compare to other materials?

	NR Latex	Nitrile	Vinyl	Copolymer
Barrier Protection	•	•	✓	✓
Fit	•	✓	✓	✓
Comfort	•	✓	✓	✓
Dexterity	•	✓	✓	✓
Tear Resistance	•	-	✓	✓
Allergy Type I	✓	•	•	•

✓ Superior • Equal - Inferior

Areas of use

- General wards.
- Labs.
- Sterilising departments.
- Pharmacy.

Glove colour

- Green coloured neoprene for easy identification as a synthetic glove.

Special features

- Delivers the protection, comfort, stretchability, fit, and tactile sensitivity of NRL.
- Powder free to eliminate powder related complications.
- Textured fingertips for better grip in wet and dry situations.
- Ambidextrous.
- Thiurams and mercaptobenzothiazoles free – reduces the risk of contact dermatitis.
- Layer pack for easier dispensing and reduction of waste – save on glove wastage.
- Size range: XS-XL (five sizes available).

Surface finish

- Textured surface for improved wet and dry grip.

Glove type

- Polymer coated powder free glove.

Cuff features

- Beaded cuff.

Glove properties

- Material: 100 per cent neoprene.
- Colour: Green.
- Donning agent: Polymer coating.
- Design: Ambidextrous with beaded cuff.
- Fingertip thickness (single): 0.16mm.
- Cuff thickness (single): 0.11mm.
- Length: 240mm.

Quality standards

- Statistically tested to an AQL of 1.5 for pinholes (exceeds world standards).
- Meets or exceeds the specifications in ASTM D 3578, Type II.
- Manufactured within the quality guidelines of ISO 9001 and CE certification.

Ordering information

Product	Code	Size	Packaging
Micro-Touch Affinity	34012001	Extra Small	100 pieces per box/ 10 boxes per carton
	34012002	Small	
	34012003	Medium	
	34012004	Large	
	34012006	Extra Large	



Ansell Limited

Level 3, 678 Victoria Street
Richmond, VIC 3121 Australia

Telephone 61 3 9270 7270
Facsimile 61 3 9270 7300

protection@ap.ansell.com
www.professional.ansell.com.au

Ansell and Micro-Touch® Affinity® are trademarks owned by Ansell Limited, or one of its affiliates.
© 2007 All Rights Reserved.

Ansell