



## COVID-19 OR FLU?

NEW SELF TEST ALLOWS PATIENTS  
TO TEST IN THE COMFORT OF HOME



### PANBIO™ COVID-19/FLU A&B PANEL SELF TEST



For patients presenting symptoms  
or suspected of current **COVID-19, Flu  
A and/or B infection**



Uses single, patient-friendly nasal swab



Suitable for all ages 1 and above\*



Results in 15 minutes



Meets more stringent CE-Mark  
IVDR requirements




Reduced environmental impact with up  
to 33% smaller and 20% lighter kit boxes  
that make shipping and storage easier\*\*

**ABBOTT HAS A RICH LEGACY AS A TRUSTED BRAND IN MORE THAN 100 COUNTRIES.**

PROUDLY DISTRIBUTED BY **EBOS HEALTHCARE**, 0800 105 501, [eboshealthcare.co.nz](https://eboshealthcare.co.nz)




## DESIGNED FOR EASY TESTING

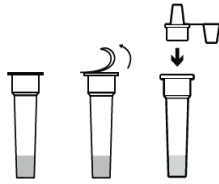


- Pre-filled buffer tube with integrated dropper cap. No manual filling needed

Nasal swab

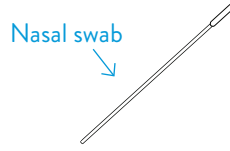


- Only one nasal swab specimen needed

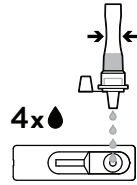

- Simple test cassette design with one sample well

## SIMPLIFIED TEST PROCEDURE

Consult Instructions for Use for complete procedure.

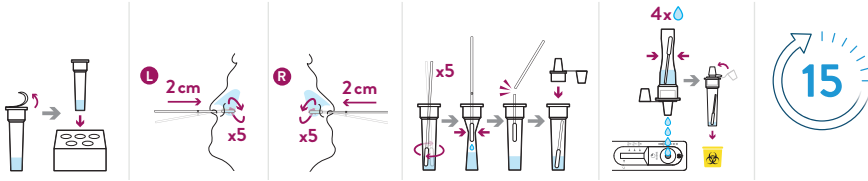

## SPECIFICATIONS

- **TEST TIME:** 15 MINUTES
- **STORAGE:** 2–30°C (shelf life 24 months)
- **CE MARK:** complies with new, more stringent CE Mark IVDR requirements
- **SAMPLE TYPE:** Nasal swab sample

## HOW ACCURATE IS THE PANBIO COVID/FLU A&B PANEL SELF TEST?

The **Panbio™ COVID-19/Flu A&B Panel Self Test** is a high-performing lateral flow test with demonstrated sensitivity for COVID-19 of 94.4% in Ct values ≤ 24 and overall specificity of 99.7% for COVID-19. For Influenza, the **Panbio™ Covid-19/Flu A&B Panel Self Test** showed 96.3% accuracy for detecting Flu A and Flu B and overall specificity of 98.6% for Flu A and Flu B.

## SUSTAINABILITY

### PANBIO COVID-19 AG SELF-TEST – OLD PACKAGING

- 220x68x30mm



### PANBIO COVID-19/FLU A&B PANEL SELF TEST – NEW PACKAGING

- 120x55x25mm



The included Nasal Swab is now 7cm shorter, allowing for smaller kit sizes and reduced waste.

15cm



8cm

Special “break-off” swab stays contained in tube, minimizing exposure.



**INTENDED PURPOSE:** Panbio™ COVID-19/Flu A&B Panel Self Test is a manual test and is intended to be used by untrained lay users (self testing) in a home setting to aid in the diagnosis of an active COVID-19, Flu A or Flu B infection.

## ORDER INFORMATION

- **Panbio™ COVID-19/Flu A&B Panel Self Test**
- **CATALOG NUMBER:**  
(1T/4T): 62FK51/62FK71

## CONTACT YOUR LOCAL REPRESENTATIVE TODAY.

[globalpointofcare.abbott](https://globalpointofcare.abbott)

\* Individuals below 16 years of age and person with reduced physical, sensory or mental capabilities should be swabbed and be supervised by an appropriate person responsible for their safety.

\*\*Compared to the Panbio™ COVID-19/Flu A&B Rapid Panel (Nasopharyngeal).

© 2024 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners. Any photos displayed are for illustrative purposes only. COL-25075-01 06/24

PROUDLY DISTRIBUTED BY **EBOS HEALTHCARE**, 0800 105 501, [eboshealthcare.co.nz](https://eboshealthcare.co.nz)

