


 **smith&nephew**  
**LEUKOSTRIP<sup>®</sup>**  
Adhesive Hypoallergenic  
Wound Closure Strips

Minimise scars, risk of allergic  
reaction and pain on removal

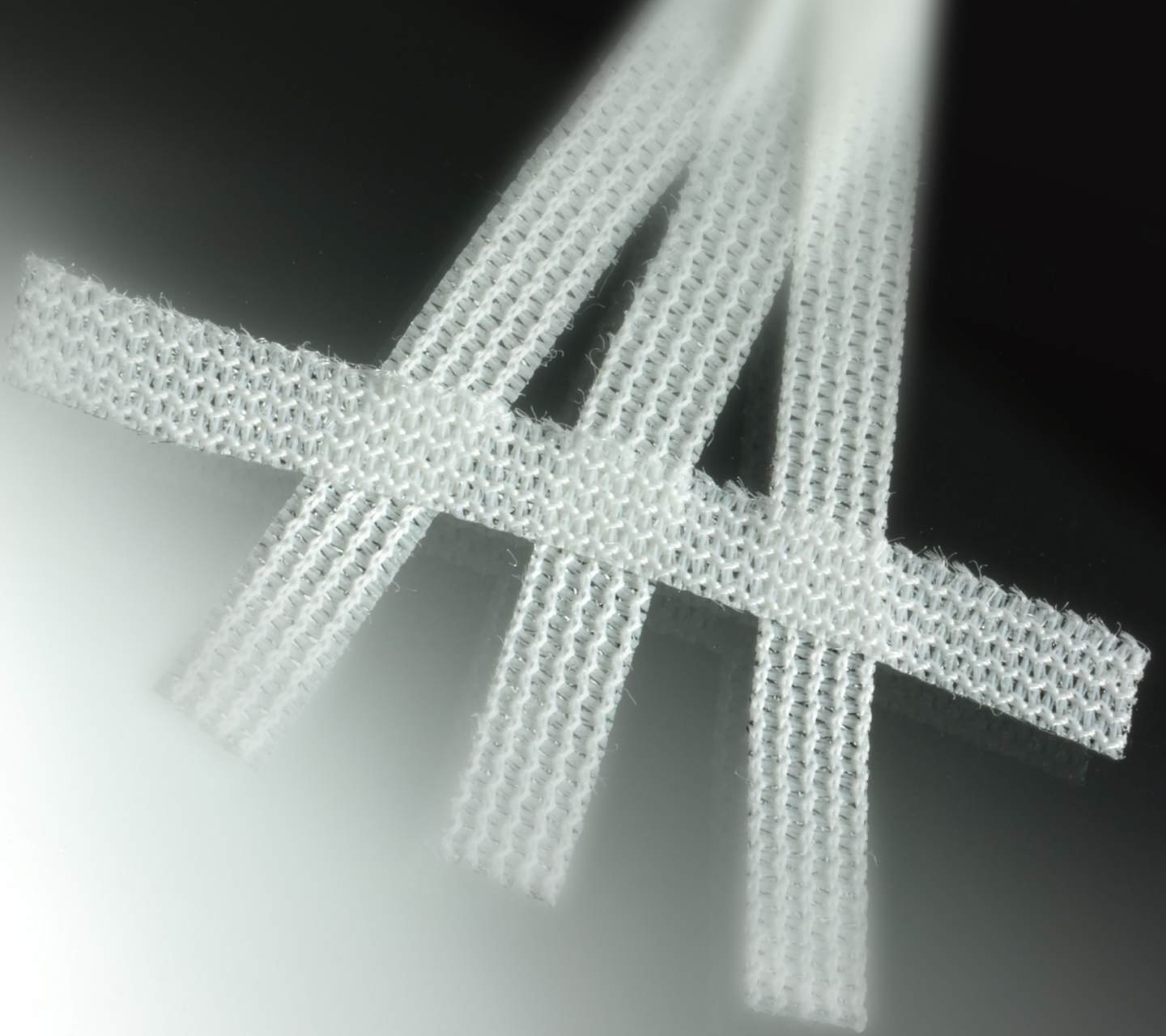

# LEUKOSTRIP<sup>◇</sup> - flexible, adhesive, hypoallergenic wound closure strips

LEUKOSTRIP wound closure strips are made of 100% polyamide material coated with a hypoallergenic adhesive. By replacing skin clips with LEUKOSTRIP on the third post-operative day, scar formation at the puncture sites is minimised<sup>1</sup>. LEUKOSTRIP can also be combined with subcutaneous or intracutaneous sutures as well as being used to supplement clip or individual over-and-over sutures.

## Advantages

### Helps protect the wound against tensile and shear forces

LEUKOSTRIP does not traumatise the edges of the wound and allows an almost pain free closure of the wound as well as minimising scar formation.

### Hypoallergenic

The hypoallergenic adhesive reduces the risk of allergic reactions.

### Reliable adhesive strength

LEUKOSTRIP adheres reliably and securely.

### Permeable

LEUKOSTRIP is permeable to water vapour and air and thus reduces the risk of maceration.

### Flexible

The elastic polyamide material returns to its original length after closing the wound.

## Indications

Primary and secondary closure of wounds.  
Replacing skin sutures following subcutaneous suturing.  
Support and relief for intracutaneous sutures, clip stitching or individual over-and-over sutures.

### Ordering Information

Code	Product description	Items per unit
<b>White</b>		
66002876	4mm x 38mm (8 strips/pouch)	Box/50
66002878	6.4mm x 76mm (3 strips/pouch)	Box/50
66002879	6.4mm x 102mm (5 strips/pouch)	Box/50
66002880	13mm x 102mm (6 strips/pouch)	Box/50
66002881	26mm x 102mm (4 strips/pouch)	Box/25
66002903	4mm x 76mm (4 strips/pouch)	Box/50
66002952	6.4mm x 76mm (3 strips/pouch)	Box/10*
66002953	13mm x 102mm (6 strips/pouch)	Box/10*
<b>Tan</b>		
66002882	4mm x 38mm (4 strips/pouch)	Box/50
66002883	6.4mm x 76mm (3 strips/pouch)	Box/50

\*Not available in New Zealand

### References

1. Lepsien G, et al. Wound healing and cosmetic results after early removal of staples and replacement with closure strips. *Akt. Chir.* 1992; 27: 60-64.



Smith & Nephew Pty Ltd  
Australia  
T +61 2 9857 3999 F +61 2 9857 3900  
www.smith-nephew.com/australia  
Customer Service  
T 13 13 60 F 1800 671 000

Smith & Nephew Ltd  
New Zealand  
T +64 9 820 2840 F +64 9 820 2841  
www.smith-nephew.com/new-zealand  
Customer Service  
T 0800 657 799 (Surgical) T 0800 807 663 (Wound)

Please consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.

◇Trademark of Smith & Nephew  
SNI2856 (11/16)